

Mercodia Ultrasensitive Rat Insulin ELISA

Directions for Use

10-1251-01
REAGENTS FOR 96 DETERMINATIONS





10-1251-10
REAGENTS FOR 10 X 96 DETERMINATIONS

For Research Use only
Not for Use in Diagnostic Procedure

Manufactured by

Mercodia AB, Sylveniusgatan 8A,
SE-754 50 Uppsala,
Sweden

EXPLANATION OF SYMBOLS USED ON LABELS

 $\Sigma = 96$	Reagents for 96 determinations
	Expiry date
	Store between 2–8°C
	Lot No.

INTENDED USE

Mercodia Ultrasensitive Rat Insulin ELISA provides a method for the quantitative determination of insulin in rat serum or plasma.

PRINCIPLE OF THE PROCEDURE

Mercodia Ultrasensitive Rat Insulin ELISA is a solid phase two-site enzyme immunoassay. It is based on the direct sandwich technique in which two monoclonal antibodies are directed against separate antigenic determinants on the insulin molecule. During incubation insulin in the sample reacts with peroxidase-conjugated anti-insulin antibodies and anti-insulin antibodies bound to the microtitration well. A simple washing step removes unbound enzyme labeled antibody. The bound conjugate is detected by reaction with 3,3',5,5'-tetramethylbenzidine. The reaction is stopped by adding acid to give a colorimetric endpoint that is read spectrophotometrically.

WARNINGS AND PRECAUTIONS

- For research use only. Not for use in diagnostic procedures.
- Not for internal or external use in humans or animals.
- The content of this kit and their residues must not be allowed to come into contact with ruminating animals or swine.
- The Stop Solution in this kit contains 0.5M H₂SO₄. Follow routine precautions for handling hazardous chemicals.

MATERIAL REQUIRED BUT NOT PROVIDED

- Pipettes for 25, 50, 100 and 200 µL (Repeating pipettes preferred for addition of enzyme conjugate solution, Substrate TMB and Stop Solution.)
- Microplate reader with 450 nm filter
- Wash device for microplates
- Tube (10–100 mL) for preparation of enzyme conjugate solution
- 1000 mL/10 L bottle
- Redistilled water
- Plate shaker (The recommended velocity is 700-900 cycles per minute, orbital movement)

REAGENTS FOR 1 X 96 KIT

Each Mercodia Ultrasensitive Rat Insulin ELISA kit (10-1251-01) contains reagents for 96 wells, sufficient for 42 samples and one Calibrator curve in duplicate. For larger series of assays, use pooled reagents from packages bearing identical lot numbers. The expiry date for the complete kit is stated on the outer label. The recommended storage temperature is 2–8°C.

Coated Plate Mouse monoclonal anti-insulin For unused microplate strips, reseal the bag using adhesive tape, store at 2–8°C and use within 8 weeks.	1 plate	96 wells 8-well strips	Ready for Use
Calibrators 1, 2, 3, 4, 5 Color coded yellow Concentration stated on vial label.	5 vials	1000 µL	Ready for Use
Calibrator 0 Color coded yellow	1 vial	5 mL	Ready for Use
Enzyme Conjugate 11X Peroxidase conjugated mouse monoclonal anti-insulin	1 vial	1.3 mL	Preparation, see below
Enzyme Conjugate Buffer Color coded blue.	1 vial	13 mL	Ready for use
Wash Buffer 21X Storage after dilution: 2–8°C for 4 weeks.	1 bottle	50 mL	Dilute with 1000 mL redistilled water to make wash buffer 1X solution.
Substrate TMB Colorless solution <i>Note! Light sensitive!</i>	1 bottle	22 mL	Ready for Use
Stop Solution 0.5 M H ₂ SO ₄	1 vial	7 mL	Ready for Use

Preparation of enzyme conjugate 1X solution

Prepare the needed volume of enzyme conjugate solution by dilution of Enzyme Conjugate 11X (1+10) in Enzyme Conjugate Buffer according to the table below. When preparing enzyme conjugate 1X solution for the whole plate, pour all of the Enzyme Conjugate Buffer into the Enzyme Conjugate 11X vial. Mix gently.

Number of strips	Enzyme Conjugate 11X	Enzyme Conjugate Buffer
12 strips	1 vial	1 vial
6 strips	600 µL	6.0 mL
4 strips	400 µL	4.0 mL

REAGENTS FOR 10 X 96 KIT

Each Mercodia Ultrasensitive Rat Insulin ELISA kit (10-1251-10) contains reagents for 10 x 96 wells, sufficient for 42 samples and one calibrator curve in duplicate on each plate. For larger series of assays, use pooled reagents from packages bearing identical lot numbers. The expiry date for the complete kit is stated on the outer label. The recommended storage temperature is 2–8°C.

Coated Plate Mouse monoclonal anti-insulin For unused microplate strips, reseal the bag using adhesive tape, store at 2–8°C and use within 8 weeks.	10 plate	96 wells 8-well strips	Ready for Use
Calibrators 1, 2, 3, 4, 5 Color coded yellow Concentration stated on vial label.	5 vials	1000 µL	Ready for Use
Calibrator 0 Color coded yellow	1 vial	5 mL	Ready for Use
Enzyme Conjugate 11X Peroxidase conjugated mouse monoclonal anti-insulin	1 vial	12 mL	Preparation, see below
Enzyme Conjugate Buffer Color coded blue.	1 vial	120 mL	Ready for use
Wash Buffer 21X	2 bottles	200 mL	Preparation, see below
Substrate TMB Colorless solution <i>Note! Light sensitive!</i>	1 bottle	220 mL	Ready for Use
Stop Solution 0.5 M H ₂ SO ₄	1 vial	70 mL	Ready for Use

Preparation of enzyme conjugate 1X solution

Prepare the needed volume of enzyme conjugate 1X solution by dilution of Enzyme Conjugate 11X in Enzyme Conjugate buffer 1+10 according to the table below. Mix gently.

Number of plates	Enzyme Conjugate 11X	Enzyme Conjugate Buffer
10 plates	1 vial	1 vial
5 plates	6000 µL	60 mL
3 plates	3600 µL	36 mL
2 plates	2400 µL	24 mL
1 plate	1200 µL	12 mL

Storage after dilution: 2–8°C for two months.

Preparation of wash buffer 1X solution

Prepare the needed volume of wash buffer by dilution of Wash Buffer 21X in redistilled water 1+20 according to the table below. Mix gently.

Number of plates	Wash Buffer 21X	Redistilled water
10 plates	2 bottles	8000 mL
5 plates	180 mL	3600 mL
3 plates	110 mL	2200 mL
2 plates	70 mL	1400 mL
1 plate	35 mL	700 mL

Storage after dilution: 2-8°C for 8 weeks.

SPECIMEN COLLECTION AND HANDLING

Serum

Collect blood by venipuncture, allow to clot, and separate the serum by centrifugation. Samples can be stored at 2–8°C up to 24 hours. For longer periods store samples at -20°C. Avoid repeated freezing and thawing.

Plasma

Collect blood by venipuncture into tubes containing heparin or EDTA as anticoagulant, and separate the plasma fraction. Samples can be stored at 2–8°C up to 24 hours. For longer periods store samples at -20°C. Avoid repeated freezing and thawing.

Preparation of samples

No dilution is normally required, however, samples containing >1.0 µg/L should be diluted 1/10 v/v with Calibrator 0. *Note!* Buffers containing sodium azide (NaN₃) can not be used for sample dilution.

TEST PROCEDURE

Prepare a calibrator curve for each assay run. All reagents and samples must be brought to room temperature before use.

1. Prepare enzyme conjugate 1X solution (according to the table on previous pages), wash buffer 1X solution and samples.
2. Prepare sufficient microplate wells to accommodate Calibrators and samples in duplicate.
3. Pipette 25 μ L each of Calibrators and samples into appropriate wells.
4. Add 100 μ L of enzyme conjugate 1X solution into each well.
5. Incubate on a plate shaker (700-900 rpm) for 2 hours at room temperature (18-25°C).
6. Wash 6 times with 700 μ L wash buffer 1X solution per well using an automatic plate washer with overflow-wash function. Do not include soak step in washing procedure.
Or manually,
Discard the reaction volume by inverting the microplate over a sink. Add 350 μ L wash solution to each well. Discard the wash solution, tap firmly several times against absorbent paper to remove excess liquid. Repeat 5 times. Avoid prolonged soaking during washing procedure.
7. Add 200 μ L Substrate TMB into each well.
8. Incubate 15 minutes at room temperature (18-25°C).
9. Add 50 μ L Stop Solution to each well.
Place the plate on the shaker for approximately 5 seconds to ensure mixing.
10. Read optical density at 450 nm and calculate results.
Read within 30 minutes.

Note! To prevent contamination between the conjugate and substrate, separate pipettes are recommended.

INTERNAL QUALITY CONTROL

Commercial controls such as Mercodia Diabetes-Antigen Control Rat, Mouse (L, M, H) (10-1220-01) and/or internal serum pools with low, intermediate and high insulin concentrations should routinely be assayed as samples, and results charted from day to day. It is good laboratory practice to record the following data for each assay: kit lot number, dilution and/or reconstitution dates of kit components, OD values for the blank, Calibrators and controls.

CALCULATION OF RESULTS

Computerized calculation

The concentration of insulin is obtained by computerized data reduction of the absorbance for the Calibrators, except for Calibrator 0, versus the concentration using cubic spline regression.

Manual Calculation

1. Plot the absorbance values obtained for the Calibrators, except for Calibrator 0, against the insulin concentration on a lin-lin paper and construct a calibrator curve.
2. Read the concentration of the samples from the calibrator curve.


Example of results

Wells	Identity	A ₄₅₀	Mean conc. µg/L
1A–B	Calibrator 0	0.077/0.070	
1C–D	Calibrator 1*	0.087/0.083	
1E–F	Calibrator 2*	0.119/0.115	
1G–H	Calibrator 3*	0.286/0.276	
2A–B	Calibrator 4*	1.063/1.027	
2C–D	Calibrator 5*	2.989/3.033	
2E–F	Sample 1	0.216/0.210	0.12
2G–H	Sample 2	0.338/0.339	0.18
3A–B	Sample 3	0.376/0.369	0.19

* Concentration stated on vial label.

Calibrator curve

A typical calibrator curve is shown here. Do not use this curve to determine actual assay results.

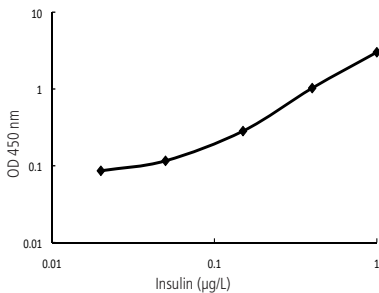

Conversion factor

1 µg corresponds to 174 pmol.

LIMITATIONS OF THE PROCEDURE

Performance limitations

Grossly lipemic, icteric or haemolysed samples do not interfere in the assay. Insulin is, however degraded over time in haemolysed samples. The degradation could give falsely low values and contribute to higher inter assay variation.

EXPECTED VALUES

Good practice dictates that each laboratory establishes its own expected range of values.

PERFORMANCE CHARACTERISTICS

Detection limit

Detection limit is defined as the Capability of Detection according to ISO11843-Part 1. Capability of Detection should be seen as a part of a method validation, rather than the lowest concentration that can be measured.

The detection limit is ≤ 0.020 µg/L as determined with the methodology described in ISO11843-Part 4.

Concentration of samples with absorbance below Calibrator 1 should not be calculated, instead expressed less or equal to (\leq) the concentration indicated on the vial for Calibrator 1.

Hook effect

Samples with a concentration up to at least 450 µg/L can be measured without giving falsely low results.

Precision

Each sample was analyzed in four replicates on 22 different occasions.

Sample	Mean value µg/L	Coefficient of variation		
		within assay %	between assay %	total assay %
1	0.11	2.9	4.8	5.0
2	0.17	2.0	4.2	4.3
3	0.19	2.6	3.5	3.8

Specificity

Human insulin	167 %
Human proinsulin	75 %
Human C-peptide	< 0.05 %
Insulin lispro (Humalog®)	167 %
IGF-I	< 0.02 %
IGF-II	< 0.02 %
Rat proinsulin I	8 %
Rat proinsulin II	51 %
Mouse proinsulin I	33 %
Mouse proinsulin II	51 %
Mouse c-peptide I	<0.002 %
Mouse c-peptide II	<0.001 %
Rat c-peptide I	<0.03 %
Rat c-peptide II	<0.03 %
Mouse insulin	75 %
Porcine insulin	476 %
Sheep insulin	179 %
Bovine insulin	78 %

WARRANTY

The performance data presented here was obtained using the procedure indicated. Any change or modification in the procedure not recommended by Mercodia AB may affect the results, in which event Mercodia AB disclaims all warranties expressed, implied or statutory, including the implied warranty of merchantability and fitness for use.

Mercodia AB and its authorised distributors, in such event, shall not be liable for damages indirect or consequential.

REFERENCES

- Korner J, Savontaus E, Chua SC, Jr., Leibel RL, Wardlaw SL (2001) Leptin regulation of AgRP and Npy mRNA in the rat hypothalamus. *J Neuroendocrinol* 13:959-966
- Kullin M, Li Z, Bondo Hansen J, Welsh N, Karlsson FA, Sandler S (2003) Protection of rat pancreatic islets by potassium channel openers against alloxan, sodium nitroprusside and interleukin-1beta mediated suppression--possible involvement of the mitochondrial membrane potential. *Diabetologia* 46:80-88
- Olsson R, Carlsson PO (2005) Better vascular engraftment and function in pancreatic islets transplanted without prior culture. *Diabetologia* 48:469-476
- Rakatzi I, Mueller H, Ritzler O, Tennagels N, Eckel J (2004) Adiponectin counteracts cytokine- and fatty acid-induced apoptosis in the pancreatic beta-cell line INS-1. *Diabetologia* 47:249-258
- Rydgren T, Sandler S (2002) Efficacy of 1400 W, a novel inhibitor of inducible nitric oxide synthase, in preventing interleukin-1beta-induced suppression of pancreatic islet function in vitro and multiple low-dose streptozotocin-induced diabetes in vivo. *Eur J Endocrinol* 147:543-551



SUMMARY PROTOCOL SHEET
Merckodia Ultrasensitive Rat Insulin ELISA

X-O Graf/Tycker/AB

Add Calibrators, Controls and Samples	25 μ L
Add enzyme conjugate 1X solution to all wells	100 μ L
Incubate	2 hours at 18-25°C on a plate shaker
Wash plate with wash buffer 1X solution	6 times
Add Substrate TMB	200 μ L
Incubate	15 minutes
Add Stop Solution	50 μ L Shake for 5 seconds to ensure mixing
Measure A ₄₅₀	Evaluate results