

# Mercodia Ultrasensitive Mouse Insulin ELISA

Directions for Use

*Please note that lot-specific  
Calibrator concentration is  
stated on vial label*





**10-1249-01**  
**REAGENTS FOR 96 DETERMINATIONS**

For Research Use only

Manufactured by

**Mercodia AB**  
Sylveniusgatan 8A  
SE-754 50 Uppsala  
Sweden

## EXPLANATION OF SYMBOLS USED ON LABELS

 $\Sigma = 96$	Reagents for 96 determinations
	Expiry date
	Store between 2–8°C
	Lot No.

## INTENDED USE

Mercodia Ultrasensitive Mouse Insulin ELISA provides a method for the quantitative determination of insulin in mouse serum or plasma.

## PRINCIPLE OF THE PROCEDURE

Mercodia Ultrasensitive Mouse Insulin ELISA is a solid phase two-site enzyme immunoassay. It is based on the direct sandwich technique in which two monoclonal antibodies are directed against separate antigenic determinants on the insulin molecule. During incubation insulin in the sample reacts with peroxidase-conjugated anti-insulin antibodies and anti-insulin antibodies bound to the microtitration well. A simple washing step removes unbound enzyme labeled antibody. The bound conjugate is detected by reaction with 3,3',5,5'-tetramethylbenzidine. The reaction is stopped by adding acid to give a colorimetric endpoint that is read spectrophotometrically.

## WARNINGS AND PRECAUTIONS

- Not for internal or external use in humans or animals.
- The content of this kit and their residues must not be allowed to come into contact with ruminating animal or swine.
- All samples should be handled as if capable of transmitting infections.
- Each well can only be used once.
- The Stop Solution contains <5% Sulphuric acid.  
The Stop Solution is labeled:



### Danger

H318 – Causes serious eye damage.

H315 – Causes skin irritation.

P280 – Wear protective gloves. Wear eye or face protection.

P264 – Wash hands thoroughly after handling.

P302 + P352 + P362 + P364 – IF ON SKIN: Wash with plenty of soap and water. Take off contaminated clothing and wash it before reuse.

P332 + P313 - If skin irritation occurs: Get medical attention.

P305 + P351 + P338 + P310 – IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a POISON CENTER or physician.

- The Enzyme Conjugate Buffer, Cal 0, 1, 2, 3, 4, 5 and Wash Buffer contain <0.06% 5-chloro-2-methyl-4-isothiazolin-3-one and 2-methyl-2H-isothiazol-3-one (3:1). The Enzyme Conjugate Buffer, the Calibrators and Wash Buffer are labeled:



### **Warning**

H317 – May cause an allergic skin reaction.

P280 – Wear protective gloves.

P261 – Avoid breathing vapour.

P272 – Contaminated work clothing should not be allowed out of the workplace.

P302 + P352 – IF ON SKIN: Wash with plenty of soap and water.

P333 + P313 – If skin irritation or rash occurs: Get medical attention.

P501 – Dispose of contents and container in accordance with all local, regional, national and international regulations.

### **MATERIAL REQUIRED BUT NOT PROVIDED**

- Pipettes with appropriate volumes (repeating pipettes preferred for addition of enzyme conjugate solution 1X, Substrate TMB and Stop Solution)
- Tubes, beakers and cylinders for reagent preparation
- Redistilled water
- Magnetic stirrer
- Vortex mixer
- Microplate reader with 450 nm filter
- Microplate shaker (700–900 cycles per minute, orbital movement)
- Microplate washing device with overflow function (recommended but not required)

## REAGENTS

Each Mercodia Ultrasensitive Mouse Insulin ELISA kit (10-1249-01) contains reagents for 96 wells, sufficient for 42 samples and one Calibrator curve in duplicate. For larger series of assays, use pooled reagents from packages bearing identical lot numbers. The expiry date for the complete kit is stated on the outer label. The recommended storage temperature is 2–8°C.

<b>Coated Plate</b> Mouse monoclonal anti-insulin For unused microplate strips, reseal the bag using adhesive tape, store at 2–8°C and use within 2 months.	1 plate	96 wells 8-well strips	Ready for use
<b>Calibrators 1, 2, 3, 4, 5</b> Mouse insulin Color coded yellow Concentration stated on vial label	5 vials	1000 µL	Ready for use
<b>Calibrator 0</b> Color coded yellow	1 vial	5 mL	Ready for use
<b>Enzyme Conjugate 11X</b> Peroxidase conjugated mouse monoclonal anti-insulin	1 vial	1.3 mL	Preparation, see below
<b>Enzyme Conjugate Buffer</b> Color coded blue	1 vial	13 mL	Ready for use
<b>Wash Buffer 21X</b> Storage after dilution: 2–8°C for 2 months	1 bottle	50 mL	Dilute with 1000 mL redistilled water to make wash buffer 1X solution.
<b>Substrate TMB</b> Colorless solution <i>Note! Light sensitive!</i>	1 bottle	22 mL	Ready for use
<b>Stop Solution</b> 0.5 M H <sub>2</sub> SO <sub>4</sub>	1 vial	7 mL	Ready for use

## Preparation of enzyme conjugate 1X solution

Prepare the needed volume of enzyme conjugate 1X solution by dilution of Enzyme Conjugate 11X (1+10) in Enzyme Conjugate Buffer according to the table below. When preparing enzyme conjugate 1X solution for the whole plate, pour all of the Enzyme Conjugate Buffer into the Enzyme Conjugate 11X vial. Mix gently.

Number of strips	Enzyme Conjugate 11X	Enzyme Conjugate Buffer
12 strips	1 vial	1 vial
8 strips	800 $\mu$ L	8 mL
4 strips	400 $\mu$ L	4 mL

Storage after dilution: 2–8°C for 2 months.

## SPECIMEN COLLECTION AND HANDLING

### Serum

Collect blood by venipuncture, allow to clot, and separate the serum by centrifugation. Samples can be stored at 2–8°C up to 24 hours. For longer periods store samples at -20°C. Avoid repeated freezing and thawing.

### Plasma

Collect blood by venipuncture into tubes containing heparin or EDTA as anticoagulant, and separate the plasma fraction by centrifugation. Samples can be stored at 2–8°C up to 24 hours. For longer periods store samples at -20°C. Avoid repeated freezing and thawing.

## Preparation of samples

No dilution is normally required for serum and plasma samples, however, samples with a concentration above Calibrator 5 should be diluted in Calibrator 0. *Note!* Buffers containing sodium azide ( $\text{NaN}_3$ ) can not be used for sample dilution.

## TEST PROCEDURE

All reagents and samples must be brought to room temperature before use. Prepare a calibrator curve for each assay run. The product has been optimized and validated without plate sealer.

1. Prepare enzyme conjugate 1X solution (according to the table on previous page), wash buffer 1X solution and samples.
2. Prepare sufficient microplate wells to accommodate Calibrators, controls and samples in duplicate.
3. Pipette 25  $\mu\text{L}$  each of Calibrators, controls and samples into appropriate wells.
4. Add 100  $\mu\text{L}$  of enzyme conjugate 1X solution into each well.
5. Incubate on a plate shaker (700-900 rpm) for 2 hours at room temperature (18-25°C).
6. Wash 6 times with 700  $\mu\text{L}$  wash buffer 1X solution per well using an automatic plate washer with overflow-wash function, after final wash, invert and tap the plate firmly against absorbent paper. Do not include soak step in washing procedure.

Or manually,

Discard the reaction volume by inverting the microplate over a sink. Add 350  $\mu\text{L}$  wash buffer 1X solution to each well. Discard the wash solution, tap firmly several times against absorbent paper to remove excess liquid. Repeat 5 times. Avoid prolonged soaking during washing procedure.

7. Add 200  $\mu\text{L}$  Substrate TMB into each well.
8. Incubate 15 minutes on the bench at room temperature (18-25°C).
9. Add 50  $\mu\text{L}$  Stop Solution to each well.  
Place the plate on the shaker for approximately 5 seconds to ensure mixing.
10. Read optical density at 450 nm and calculate results.  
Read within 30 minutes.

*Note!* To prevent contamination between the conjugate and substrate, separate pipettes are recommended.

## INTERNAL QUALITY CONTROL

Commercial controls such as Mercodia Diabetes antigen Control, Rat and Mouse, Low, Medium & High (10-1220-01) and/or internal serum pools with low, intermediate and high insulin concentrations should routinely be assayed as samples, and results charted from day to day.

It is good laboratory practice to record the following data for each assay: kit lot number, preparation dates of kit components, OD values for the blank, Calibrators and Controls.

## CALCULATION OF RESULTS

### Computerized calculation

The concentration of insulin is obtained by computerized data reduction of the absorbance for the Calibrators, except for Calibrator 0, versus the concentration using cubic spline regression.

## Manual Calculation

1. Plot the absorbance values obtained for the Calibrators, except for Calibrator 0, against the insulin concentration on a log-log paper and construct a calibrator curve.
2. Read the concentration of the samples from the calibrator curve.


## Example of results

Wells	Identity	A <sub>450</sub>	Mean conc. µg/L
1A–B	Calibrator 0	0.099/0.094	
1C–D	Calibrator 1*	0.115/0.109	
1E–F	Calibrator 2*	0.136/0.136	
1G–H	Calibrator 3*	0.248/0.261	
2A–B	Calibrator 4*	0.703/0.730	
2C–D	Calibrator 5*	2.877/2.881	
2E–F	Sample 1	0.184/0.180	0.131
2G–H	Sample 2	0.248/0.260	0.200
3A–B	Sample 3	0.386/0.409	0.308

\* Concentration stated on vial label.

## Calibrator curve

A typical calibrator curve is shown here. Do not use this curve to determine actual assay results.

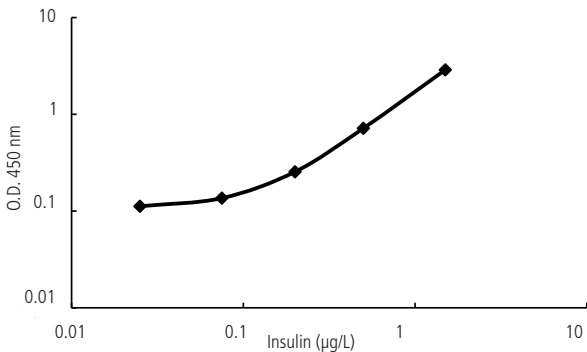



## Conversion factor

1 µg corresponds to 174 pmol

## LIMITATIONS OF THE PROCEDURE

### Performance limitations

Grossly lipemic, icteric or hemolyzed samples do not interfere in the assay. However, hemolysis in serum and plasma samples may result in a degradation of insulin which could give falsely low values and contributes to higher inter assay variation. The degradation is dependent on time, temperature and the hemoglobin concentration. Keep hemolyzed samples cold or on ice to prevent the insulin degradation.

Separate pipette tips should be used when pipetting the conjugate and the substrate.

### EXPECTED VALUES

Good practice dictates that each laboratory establishes its own expected range of values.

## PERFORMANCE CHARACTERISTICS

### Detection limit

Capability of Detection should be seen as part of a method validation, rather than the lowest concentration that can be measured. The detection limit is  $\leq 0.025$  µg/L as determined by the methodology described in ISO11843- Part 4. Concentration of samples with absorbance below Calibrator 1 should not be calculated, instead expressed as less or equal to ( $\leq$ ) the concentration indicated on the vial for Calibrator 1.

### Recovery

Recovery upon addition is 76-94 (mean 85)%.

Recovery upon dilution is 110-132 (mean 123)%.

### Hook effect

Samples with a concentration up to at least 450 µg/L can be measured without giving falsely low results.

### Precision

Each sample was analyzed in 4 replicates on 16 different occasions.

---

Sample	Mean value µg/L	Coefficient of variation		
		within assay %	between assay %	total assay %
1	0.14	2.9	4.9	5.1
2	0.22	1.5	3.9	3.9
3	0.26	2.2	4.5	4.6

---

## Specificity

Human insulin	195%
Human proinsulin	82%
Human C-peptide	n.d.
IGF-I	n.d.
IGF-II	n.d.
Rat insulin	146%
Rat proinsulin	14%
Rat C-peptide	n.d.
Porcine insulin	628%
Ovine insulin	256%
Bovine insulin	110%

n.d.= not detected

## WARRANTY

The performance data presented here was obtained using the procedure indicated. Any change or modification in the procedure not recommended by Mercodia AB may affect the results, in which event Mercodia AB disclaims all warranties expressed, implied or statutory, including the implied warranty of merchantability and fitness for use.

Mercodia AB and its authorised distributors, in such event, shall not be liable for damages indirect or consequential.

## REFERENCES

Blyszczuk P, Czyz J, Kania G, Wagner M, Roll U, St-Onge L and Wobus AM (2003) Expression of Pax4 in embryonic stem cells promotes differentiation of nestin-positive progenitor and insulin-producing cells. *Proc Natl Acad Sci U S A* 100:998-1003.

Burcelin R, Crivelli V, Perrin C, Da Costa A, Mu J, Kahn BB, Birnbaum MJ, Kahn CR, Vollenweider P and Thorens B (2003) GLUT4, AMP kinase, but not the insulin receptor, are required for hepatoportal glucose sensor-stimulated muscle glucose utilization. *J Clin Invest* 111:1555-1562.

Friesen NT, Buchau AS, Schott-Ohly P, Lgssiar A and Gleichmann H (2004) Generation of hydrogen peroxide and failure of antioxidative responses in pancreatic islets of male C57BL/6 mice are associated with diabetes induced by multiple low doses of streptozotocin. *Diabetologia* 47:676-685.

Jaeckel E, Lipes MA and von Boehmer H (2004) Recessive tolerance to preproinsulin 2 reduces but does not abolish type 1 diabetes. *Nat Immunol* 5:1028-1035.

Further references can be found on our website: [www.mercodia.com](http://www.mercodia.com)











	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C												
D												
E												
F												
G												
H												

Experiment:

Date:

KIT LOT#:

## SUMMARY PROTOCOL SHEET

### Merckodia Ultrasensitive Mouse Insulin ELISA

Add Calibrators, controls* and Samples	25 $\mu$ L
Add enzyme conjugate 1X solution to all wells	100 $\mu$ L
Incubate	2 hours at 18-25°C on a plate shaker (700-900 rpm)
Wash plate with wash buffer 1X solution	700 $\mu$ L, 6 times
Add Substrate TMB	200 $\mu$ L
Incubate	15 minutes at 18-25 °C
Add Stop Solution	50 $\mu$ L Shake for 5 seconds to ensure mixing
Measure A <sub>450</sub>	Evaluate results

\*not included

For full details see page 6

For technical support please contact: [support@merckodia.com](mailto:support@merckodia.com)